




## Extensive choice to suit your experimental needs

Many cryostats are used in experimental apparatus where the samples must be irradiated or measurements made on emitted radiation from such samples. It is, therefore, essential that appropriate windows can be incorporated into your cryostat to permit radiation to pass through the sample space.

## Window materials

The type of window material selected will be determined by the wavelength and intensity of the radiation, the environmental temperature and whether or not beam polarisation is required.

Oxford Instruments Superconductivity can offer a wide range of windows options:

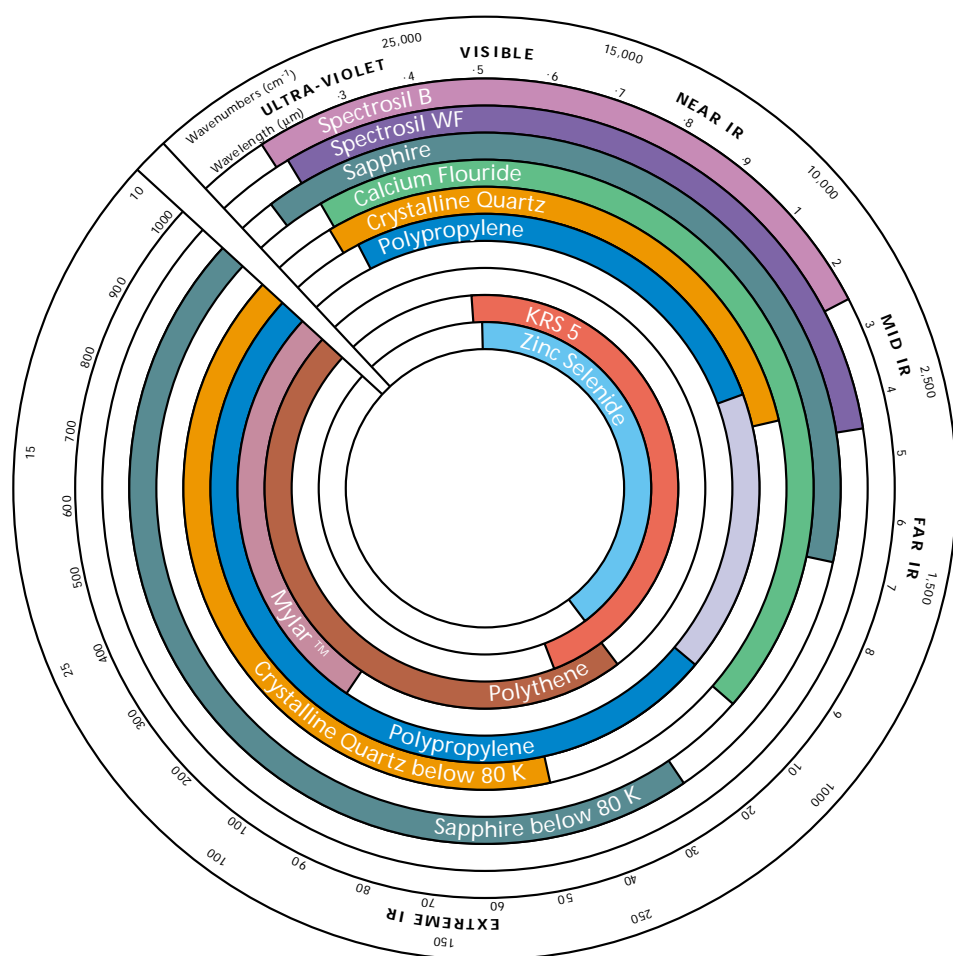


## Window mounting

Room temperature windows which must be vacuum tight are generally sealed onto a rubber 'O' ring. The mounting method minimises strain on the windows and allows easy removal if required.

Radiation shield windows do not need to be vacuum tight but it is important that they are thermally anchored. This is achieved by retaining the window in a copper mount using a spring clip.


The windows on sample cells must be vacuum tight at low temperatures and two sealing methods, indium or copper gaskets, are used depending on whether the upper temperature is limited to 300 K or 500 K respectively.



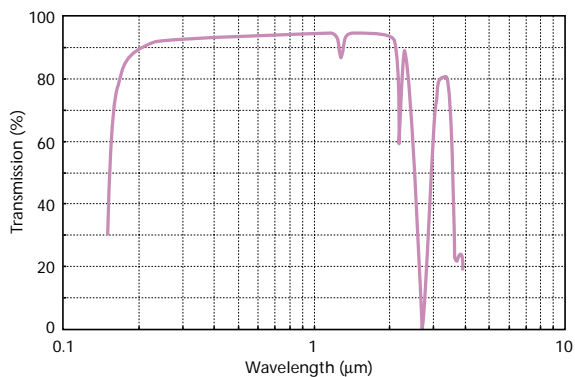


### Spectrosil B

Spectrosil B is a synthetic vitreous silica. The material itself is not birefringent and standard strain-relieved low temperature window mounts ensure that Spectrosil B windows do not affect polarised light. It is an excellent material for filtering out the near-mid IR room temperature thermal radiation, resulting in a low thermal load on the cryostat. The transmission data is for a 2 mm thick window.

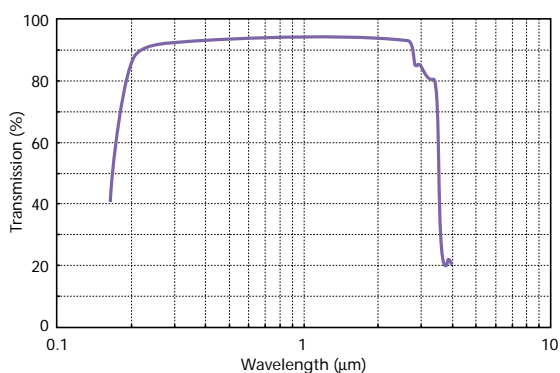


### Spectrosil WF

Spectrosil WF is a water free synthetic vitreous silica. It has the same material characteristics as Spectrosil B, but has a transmission range which extends further into the IR. The transmission data is for a 2 mm thick window.

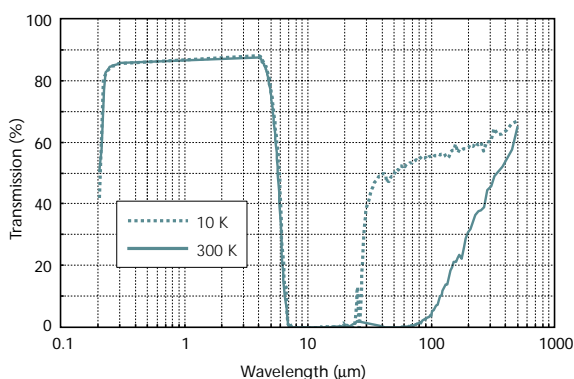


### Sapphire

Sapphire is an alternative window for visible and near IR applications. It is the only window which may be used at high temperatures (up to 500 K). At temperatures below 80 K the transmission of sapphire increases in the extreme IR region. This effect allows sapphire to be used as a cold inner window for IR measurements. Sapphire is birefringent. The transmission data is for a 2 mm thick window.

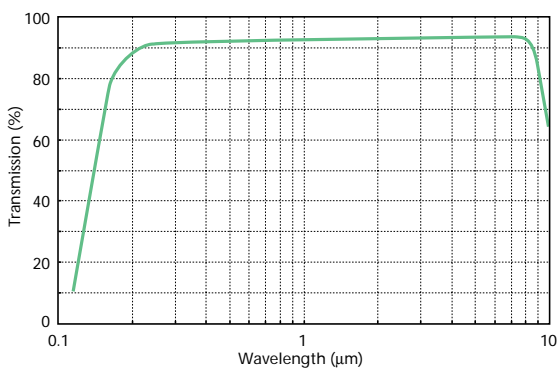


### Calcium Fluoride

Calcium fluoride may be used in the near and mid IR. It cannot be used as a cold inner window due to its mechanical and thermal properties. The transmission data is for a 2 mm thick window.

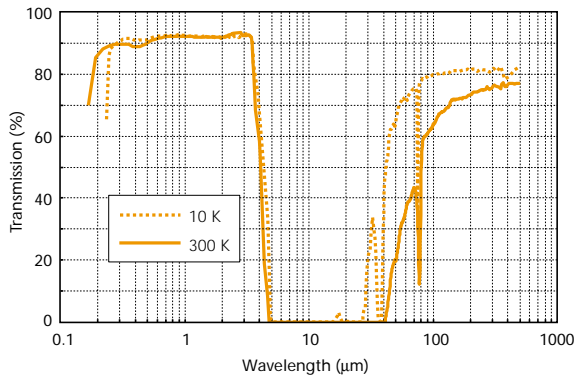


### Crystalline Quartz

Crystalline quartz is a window for the visible and near IR. It exhibits a strong temperature dependence in the far IR. At temperatures below 80 K the transmission is greater than 70% for wavelengths greater than 80 μm. This material may be used as a cold inner window for far IR applications. Crystalline quartz is birefringent and is available cut with the z-axis normal to the window to reduce polarisation effects. The transmission data is for a 2 mm thick window.

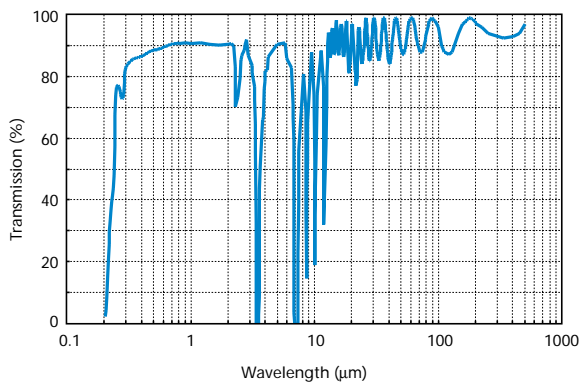


### Polypropylene

Polypropylene is a polymer with a useful transmission range which runs from 0.2 μm to the millimetre region with only a few absorption bands in the mid IR. At temperatures above 100 K this material is porous to helium gas. The transmission data is for a window 60 μm thick which may be used as a cold window to temperatures as low as 1.5 K. Polypropylene windows may be damaged if used with high power optical beams.

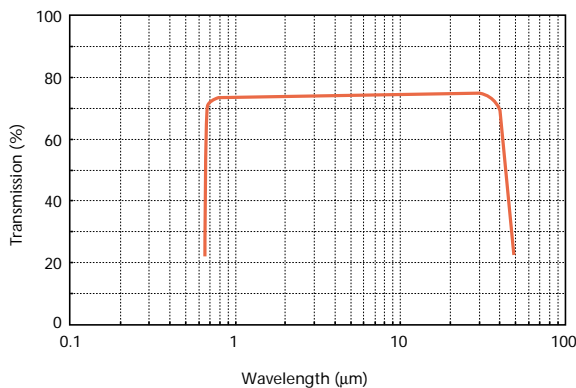


### KRS 5

KRS 5 (Thallium Bromide-Thallium Iodide) is a useful window for the mid IR. The complex nature of the cold window mount and the necessity to make it vacuum tight results in the clear access diameter being reduced by 2 mm when compared with other materials. The material is toxic and should be handled with care.

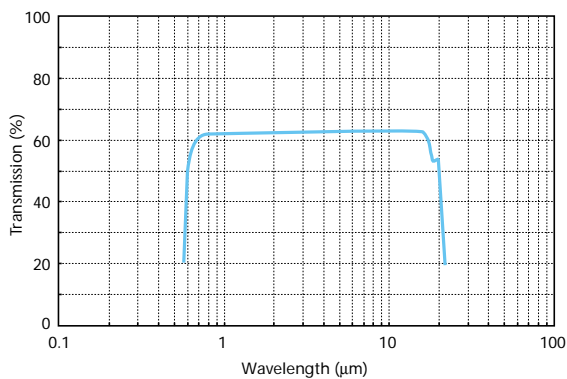


### Zinc Selenide

Zinc Selenide is a reliable cryogenic window for the mid IR. Care should be taken with handling as it is toxic and scratches easily. The transmission data is for a 2 mm thick window.

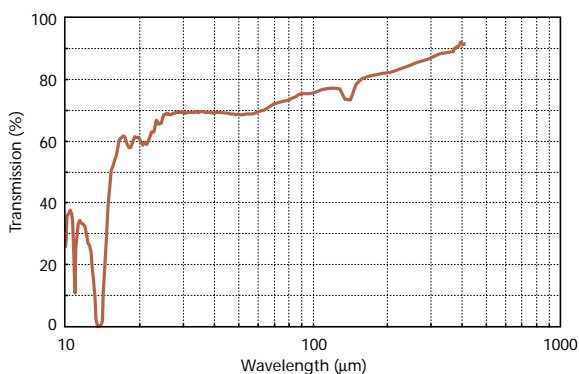

### Polythene

Polythene is a polymer which may be used in the extreme IR. The material is porous to helium gas above 100 K. The transmission data is for a 1 mm thick window.

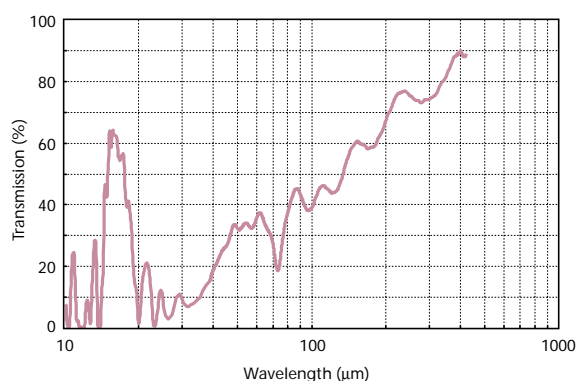

### Mylar™

Mylar™ is a polyester film. Standard windows are 125 μm thick and may be used as a cryogenic window in the extreme IR. Mylar is porous to helium gas above 100 K.

### Aluminised Mylar™

Aluminised Mylar™ is a useful window for X-ray and gamma ray applications. For X-rays its transmission is greater than 90% with some diffuse scattering at low angles. Standard windows are 125 μm thick with a 0.3 μm thick aluminium coating on one side.

Others available on request.

### Oxford Instruments Superconductivity

UK  
Tubney Woods, Abingdon  
Oxfordshire OX13 5QX  
Tel: +44 (0)1865 393 200  
Fax: +44 (0)1865 393 333

China  
Rm. 14-F, No.1 Plaza  
No.800 Nanjing East Road  
Shanghai 200001  
Tel: +86 21 63608530/1/2/3  
Fax: +86 21 63608535

Germany  
Otto-von Guericke Ring 10  
D-65205 Wiesbaden  
Tel: +49 6122 937 171  
Fax: +49 6122 937 175

Italy  
Via Leone Tolstoj 86  
20098 San Giuliano  
Milanese  
Milan  
Tel: +39 02 98 2215  
Fax: +39 02 98 2414

Japan  
Haseman Building  
2-11-6 Tomioka, Koto-ku  
Tokyo 135-0047  
Tel: +81 03 5245 3261  
Fax: +81 03 5245 4472

Spain  
Avda Mata Piñonera, 2  
28700 San Sebastian de los Reyes  
Madrid  
Tel: +34 91 659 0740  
Fax: +34 91 654 6794

U.S.A.  
130A Baker Ave. Ext.  
Concord, MA 01742-2121  
Tel: +1 978 369 9933  
Fax: +1 978 369 6616

E-mail:  
superconductivity@oxinst.co.uk

Visit our web site at  
[www.oxford-instruments.com](http://www.oxford-instruments.com)

This publication is the copyright of Oxford Instruments Superconductivity Limited and provides outline information only which (unless agreed by the company in writing) may not be used, applied or reproduced for any purpose or form part of any order or contract or be regarded as a representation relating to the products or services concerned. Oxford Instruments' policy is one of continued improvement. The company reserves the right to alter, without notice, the specification, design or conditions of supply of any product or service. Oxford Instruments acknowledges all trade marks and registrations

© Oxford Instruments Superconductivity Limited, 2003. All rights reserved.

Reference No: HCC OI634 08/03

UK sites of Oxford Instruments Superconductivity Limited, operate a Quality Management System approved to the requirements of BS EN ISO 9001.



Certificate No Q4118

